
**FEEL!<sup>2</sup>**  
**NECK**  
Shoulder Belt

---

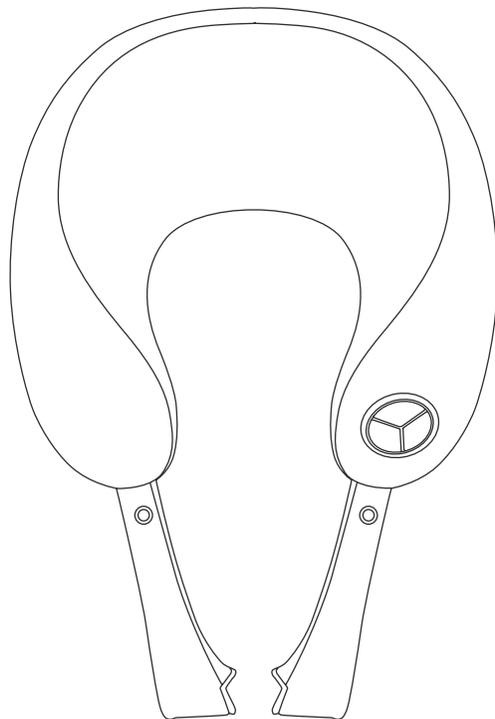100
95
75
25
5
0

## TABLE OF CONTENTS


IMPORTANT INFORMATION-----	1
PRODUCT OVERVIEW-----	2
OPERATION INSTRUCTION-----	3
TROUBLESHOOTING-----	5
PRODUCT MAINTENANCE-----	6
PRODUCT SPECIFICATION-----	7

## IMPORTANT INFORMATION

1. It is recommended to drink 250cc of water before and after using the SOFO F-Wrap II Shoulder Strap Massager for a better effect.
2. Please do not massage over 10 minutes at a single spot of the body.
3. Strictly prohibited the use of SOFO F-Wrap II Shoulder Strap Massager on an empty stomach, after meal, intoxicated on alcohol and after strenuous sports. Do not use during high fever.
4. Preliminary usage may cause slight itch to the skin, this is a normal symptom.
5. For people not sensitive to heat, please use it precautious.
6. For people with osteoporosis, heart disease, high blood pressure and under menstrual, please seek for doctor advise. Children and pregnant woman is prohibited to use.
7. For non-professional, please do not attempt to disassemble, service or repair the SOFO F-Wrap II Shoulder Strap Massager yourself.
8. If there is damage on the power cord cable, the SOFO F-Wrap II Shoulder Strap Massager must only be professionally serviced or repaired that follows the manufacturing specification for safety purpose.

## PRODUCT OVERVIEW

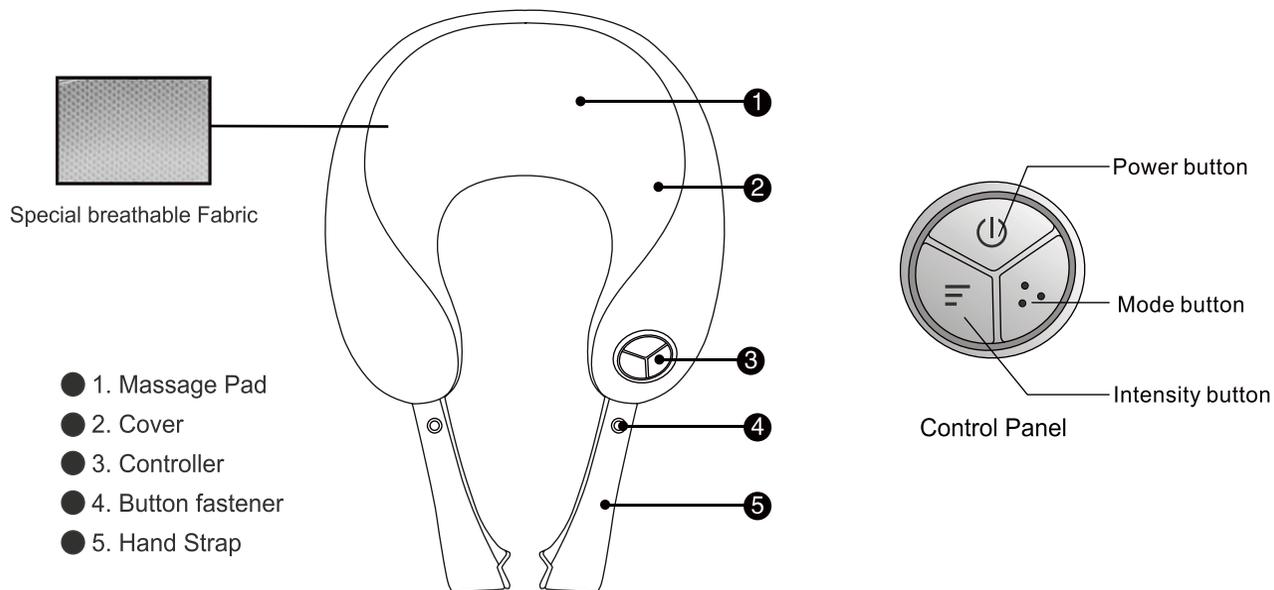 All images shown are based on the actual product and are subjected to changes without notice.


## OPERATION INSTRUCTION

1. Power button:  Power ON/OFF. Pushing the power button with a "beep" sound turns the machine ON.  
Pushing the power button again with a "beep" sound turns the machine OFF.
2. Mode button:  Mode selection control. A total of 5 modes, every push of the mode button changes the mode from one to the other. Mode 5 changes back to mode 1 when the button is pushed again and it recycles.
3. Intensity button:  Intensity selection control. A total of 5 levels, it starts with the lowest intensity (first level) when the machine is ON. every push of the intensity button changes the intensity level from one to the other. Level 5 changes back to level 1 when the button is pushed again and it recycles.
4. Every massage lasts for 10 minutes. Please push the power button and then unplug the power cord from the electrical outlet if you intended to stop the operation.

**Buzzer tips:**

The buzzer will sound upon plugging in the power and every button that is pushed.


## OPERATION MANUAL

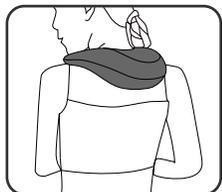
Shoulder and neck massage


Back massage

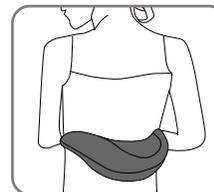
Waist massage

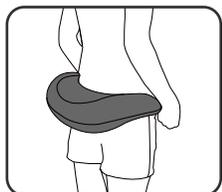Hip massage

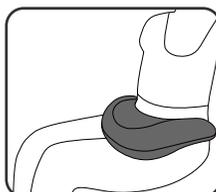Front leg massage

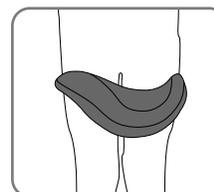Back leg massage

100
95
75
25
5
0

## TROUBLESHOOTING

MALFUNCTIONS	POSSIBLE FACTOR	SOLUTION
Product does not work after connecting to a power source.	The F-wrap II has not been turned ON?	Ensure the main switch is turned on by pushing the power button.
	Loose contact between the power plug and the electrical outlet?	Ensure the power plug is firmly connected to the electrical outlet.
Product stopped functioning halfway while using it.	Our product is set to finish or stop after 10 minutes of operation.	To continue massage, please push the power button on the controller again.
	Overuse, overheat protection.	Wait for it to cool down to start using again.
	The weight on product is overloaded, ie. somebody stepped or heavy object is loading on the massage pad.	Unload the weight on the item, remove and reconnect the AC plug to the electrical outlet.
Product produce Loud noise during using it. (motor sound)	These sounds are the result of mechanical structure.	The sound is normal.
The power cord is damaged, the product is smoking or abnormally hot.	Stop use immediately and contact for customer service or repair.	

If malfunction is still not solved with the above advice, please send it repairs or services.

## PRODUCT MAINTENANCE

1. Keep F-wrap II clean from dust and away from places of high humidity for storage.
2. Do not store in places of high temperature or places that burn easily. Do not place under direct sunlight for a long period of time.
3. Do not use cleaning agents such as thinner, oil or bleach to clean the F-wrap II.
4. Use only a soft, dry cloth to clean. Do not use wet cloth to clean the F-wrap II or the power plug cable.

## PRODUCT SPECIFICATION

Product name	Shoulder belt
Product code	FEELNECK2
Rated Voltage	AC 220V / 50Hz
Power consumption	60W
Rating time	10 minutes
Implementation of standards	GB4706.10 – 2008
to changes without notice.	

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Авторское право ©2012 г.
ERGONOVA GmbH